

## 9106 RODISIL CLOSED PATELLA KNEE BRACE WITH METAL HINGES

Closed elastic knee support with silicone pad and metal, polycentric hinges.

### Features:

Breathable elastic knee support with material made of up to 140 filaments, which provides active compression in 3 dimensions, with closed kneecap and with distal and proximal extremes designed to lighten the pressure preventing excessive compression. With finer silicone pad in the distal region to release the patellar tendon and metal, polycentric hinges on both sides to control mediolateral instabilities. **The Rodisil provides a therapeutic compression of 20.33-36 mmHg.**

### Function:

Support, compression, and adjustment. The silicone pad produces an intermittent micromassage that favours re-absorption of edemas and hematomas with every knee movement.


### Indications:

- Tendomyopathy
- First or second degree ligament sprains
- Mild Osteoarthritis
- Post-traumatic and post-operative use
- Sensation of instability (patellar and/or medial-lateral)



Bilateral design fits both left and right knee

ITEM	SIZE	THIGH CIRCUMFERENCE IN CM (A)	CALF CIRCUMFERENCE IN CM (B)
9106-XS	X-SMALL	36 - 39	26 - 29
9106-SM	SMALL	39 - 42	29 - 32
9106-MD	MEDIUM	42 - 45	32 - 35
9106-LG	LARGE	45 - 48	35 - 38
9106-XL	X-LARGE	48 - 51	38 - 41
9106-XXL	2X-LARGE	51 - 54	41 - 44
9106-XXXL	3X-LARGE	54 - 57	44 - 47
9106-4XL	4X-LARGE	57 - 60	47 - 50
9106-5XL	5X-LARGE	60 - 63	50 - 63

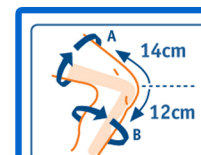

A-Point of measurement: outline halfway up thigh  
B-Point of measurement: outline halfway up calf  
Height: 34 cm