

NEW

Back n Black™


with Dual Panels



The patented design features of the Back N Black™ with Dual Panels provide your patients with the support and comfort they need in the treatment of lumbar spine related disorders.

Indications:

Chronic or Acute Low Back Pain, Degenerative Disc Disease, Lumbar Muscle Strains, Lumbar Ligamentous Sprains, Post Operative Lumbar Stabilization, Postural Support

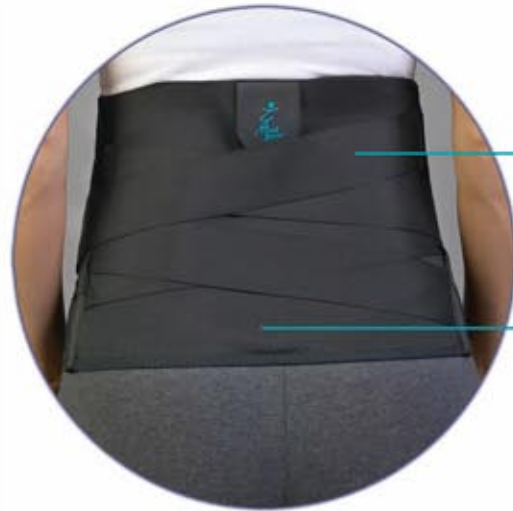

Rigid Anterior and Posterior Panels

The rigid anterior and posterior panels work in concert with the inelastic nylon strapping system to increase the intra-abdominal pressure, thereby decreasing the load on the intervertebral discs and surrounding musculoskeletal structures of the lumbar spine. Both panels may be heated and remolded to provide a custom fit.

Nylon Strapping System

The hand closure strap on the nylon strapping system allows for ease of application and tightening to create the desired level of compression.

Breathable Elastic Panels

The breathable elastic panels on the interior offer adjustability in fitting various abdominal shapes and enhance patient compliance.

Part Id	Size	Waist Circumference
163252	small	27" - 31"
163254	medium	32" - 38"
163255	large	39" - 45"
163256	x-large	46" - 52"
163257	xx-large	53" - 59"

Medical Specialties, Inc.
4600-K Lebanon Rd.
Charlotte, NC 28227
p: 800-334-4143 f: 800-694-9060
p: 704-694-2434 f: 704-694-9060
e-mail: request@medspec.com
www.medspec.com
suggested HCPCS: L0627
U.S. Patent # 5,722,940
U.S. Patent # 5,984,885

© 2006 Medical Specialties, Inc.

CE
rev 1/06