

AB01 BOXIA A.F.O. DROP FOOT



Characteristics:

The Boxia® drop-foot support raises the forefoot, allowing plantar flexion at the heel when stepping off the foot. It prevents the risk of trips over everyday items, such as rugs, curbs and irregular surfaces. The Boxia® device consists of 2 separate parts; a support above the ankles and a non-slip support section to put between the tongue and shoelaces via a central hook-type attachment (must be worn with laced shoes). Both parts are linked by an elastic strap which is adjustable according to how much bending is required. The back of the strap has interior gel padding for comfort on the Achilles tendon, while the instep area is lined to prevent unwanted friction and ease compression.

Indications:

- Stroke
- Brain Injury
- Flaccid Paralysis

How it Works:

The Boxia drop foot ankle brace's supramalleolar support stores a large amount of energy, similar to a pivot (support point), to enable it to transfer the necessary forces to lift a paralysed foot. In cases in which a patient shows intolerance due to excessive force, the Boxia calf support is recommended. It helps to stabilise the orthosis and restrict forces and stresses by transferring them from the ankle to the calf.

Benefits:

Discreet, comfortable and flexible. It occupies little space in the shoe and is easily adjusted.



ITEM	SIZE	PERIMETER IN CM
AB01-XS	X-Small / 0	13 - 17
AB01-SM	Small / 1	17 - 21
AB01-MD	Medium / 2	21 - 25
AB01-LG	Large / 3	25 - 29



See Reverse for additional accessories

ABI2 INTERIOR ARCH SUPPORT



Characteristics:

Anatomically designed to the shape of either the left or right foot, for supporting the arch area with padding.

Please note that part ABI2 contains only the lower part for the foot. You must therefore already have part AB01 (Boxia dropfoot support).

RIGHT	LEFT	SIZE	PERIMETER IN CM
AB12D-SM	AB12I-SM	Small / 1	17 - 21
AB12D-MD	AB12I-MD	Medium / 2	21 - 25
AB12D-LG	AB12I-LG	Large / 3	25 - 29





Measure
around the
instep



ABI I BOXIA ATTACHMENT HOOK

This additional Hook Attachment is one size only and is available on special order only.

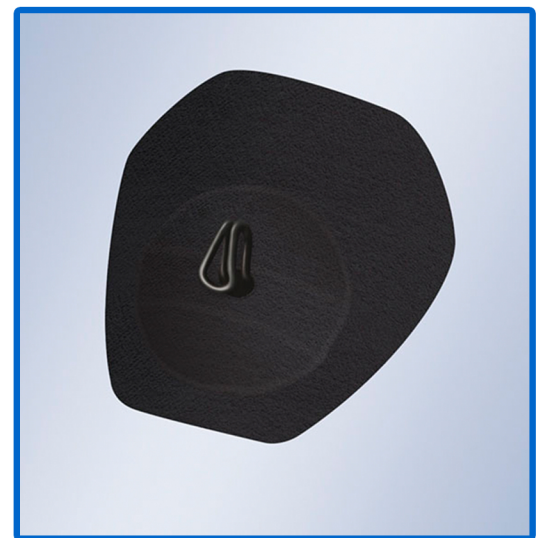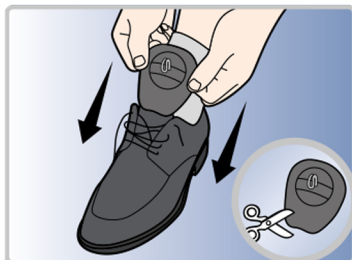

See Reverse for additional accessories