


# ASO<sup>®</sup> Ankle Stabilizing Orthosis



## The Choice of Champions<sup>SM</sup>

With an improved boot design for a better fit, the ASO ankle stabilizer maintains its leadership position in the ankle brace market. Give your athletes and patients the proven technology of the ASO ankle stabilizer. The leader continues to stay several steps ahead...

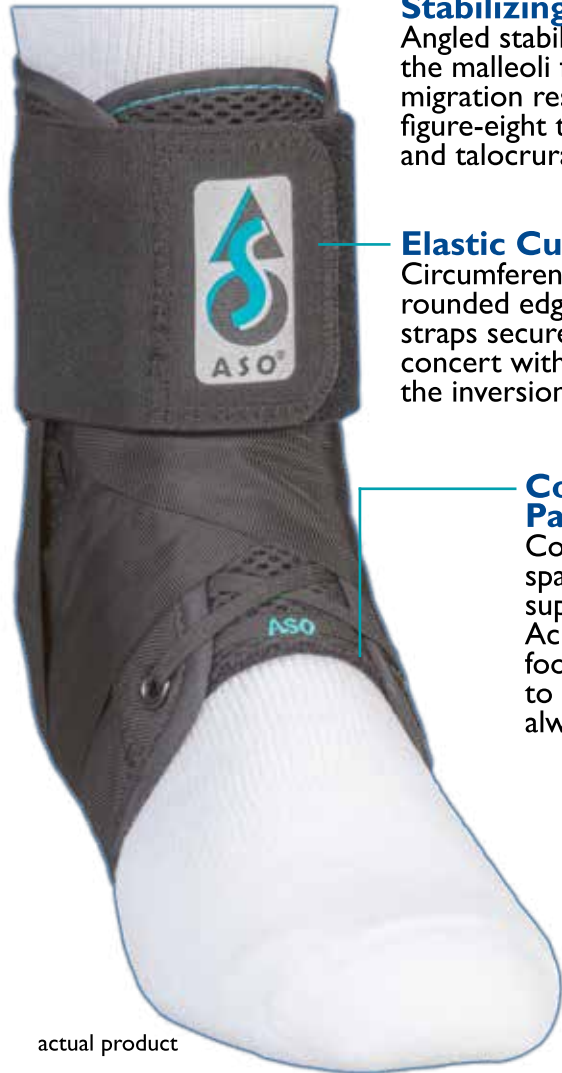

actual product

### Stabilizing Straps

Angled stabilizing straps contour above the malleoli for enhanced comfort and migration resistance. They form a complete figure-eight to help stabilize the subtalar and talocrural joints.

### Elastic Cuff

Circumferential woven elastic cuff with rounded edges keeps laces and stabilizing straps secure. In addition, it works in concert with the upper boot to resist the inversion or eversion forces.

### CoolFlex Achilles Pad & Tongue

CoolFlex<sup>™</sup> padding contains spandex fiber to provide superior comfort along the Achilles tendon and dorsum of foot. The shoe lace is attached to the tongue to ensure lace is always centered.



Elastic cuff has been removed and white straps used in the above photo for illustrative purposes.

### Invisible Seam

For improved comfort and durability, the seam is now positioned along the rear portion of the boot. There is no longer a seam under the foot.

### Indications:

- Inversion or Eversion Ankle Sprains
- Medial or Lateral Ankle Instability
- Subtalar Joint Instability

PDAC Verified HCPCS:  
L1902

Men's Shoe Size	6-7	7 1/2-9	9 1/2-11	11 1/2-13	13 1/2-15	15+
Narrow Foot	pedi	xs	s	m	l	xl
Regular Foot	xs	s	m	l	xl	xxl
Wide Foot	s	m	l	xl	xxl	xxxl
Women's Shoe Size	5-6	6 1/2-8	8 1/2-10	10 1/2-12	12 1/2-14	
Narrow Foot	pedi	xs	s	m	l	
Regular Foot	xs	s	m	m	l	
Wide Foot	s	m	l	l	xl	

Bilateral design fits left or right foot

Black	White	Size	Ankle Circum.*
264019	264009	pedi	9" - 10"
264011	264001	x-small	10" - 11"
264012	264002	small	11" - 12"
264014	264004	medium	12" - 13"
264015	264005	large	13" - 14"
264016	264006	x-large	14" - 15"
264017	264007	xx-large	15" - 16"
264018	264008	3x-large	16" +

\*measure ankle circumference according to diagram below



Medical Specialties, Inc. 4600-K Lebanon Rd. Charlotte, NC 28227  
p: 800-334-4143 f: 800-694-9060 p: 704-694-2434 f: 704-694-9060  
e-mail: request@medspec.com www.medspec.com

Patents: US 8,808,215; US 9,375,339; US 10,383,404 and other patents pending  
ASO is a registered trademark of Medical Specialties, Inc. © 2009-2023 Medical Specialties, Inc.

CE  
rev 02/23