





Professional Orthopedic Products

## CMCcare Thumb Brace

This thumb brace is suitable for daily and / or long-term use and is characterized by a high degree of efficiency, comfort and quality with the objective to maintain maximum hand function.

The CMCcare is designed for optimum support of the thumb saddle joint and positions the metacarpal of the thumb in a stable, functional position with a good gripper function. The anatomical design with rounded edges, the simple locking system and materials ensure a good fit and optimum comfort.

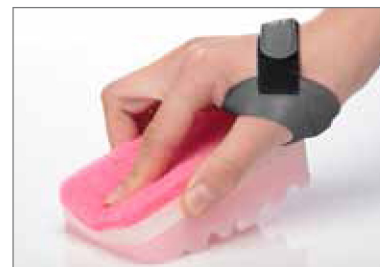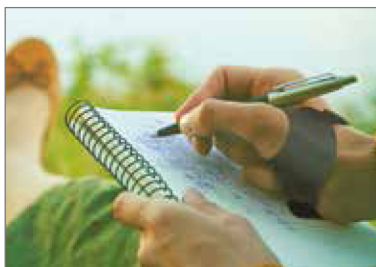

### Indications

- Osteoarthritis (thumb) (ICD 10: M19.94)
- Instability of the thumb basal joint (CMC-1)
- Postoperative after arthroplasty of the CMC-1 joint

### Product function

- Reducing range of motion joints
- Preventing deformity

### Features

- Anatomically shaped
- Integrated shapeable aluminum core
- Adjustable Velcro closure
- Additional Velcro strap
- Material: TPU plastic, aluminum, nylon. No Latex.
- Color: gray/black



### Size Chart

Item number		Unit	Size	Circumference	
Right	Left				
8712R-S	8712L-S	Piece	Small	16,0 - 19,5 cm	
8712R-M	8712L-M	Piece	Medium	19,5 - 22,5 cm	
8712R-L	8712L-L	Piece	Large	22,5 - 26,0 cm	

Measurement: to select the correct size, measure the circumference around the center hand (without thumb)

Available at: Professional Orthopedic Products  
 Tel: 1-866-334-4441  
 Fax: 1-888-484-4431  
 E-mail: [csr@poportho.com](mailto:csr@poportho.com)  
[www.poportho.com](http://www.poportho.com)